

HARMONY

Vibroacoustic Lounger



Where sound, vibration and stillness come together to restore body and mind.

From clinics to wellness spaces, the Harmony Vibroacoustic Lounger combines integrated sound and vibration to support recovery, relaxation and elevated wellness experiences.

Designed to go beyond traditional relaxation, low-frequency vibration is experienced through the body while full body cushioning helps create a deeper sense of calm without disrupting the surrounding space.





CORPORATE WELLNESS

Let employees relax while technology restores balance, reducing stress and enhancing performance in real-time. This wellness service helps them reach their full potential with less stress.



RESIDENTIAL USE

Cutting-edge technology for your family. In just 15 minutes, it helps reduce stress by up to 54%, improve sleep, ease full-body inflammation, and promote deep recovery. Vibroacoustic Therapy uses medical-grade vibrations and neural music to recalibrate the body. Perfect addition to your home Gym or general living spaces.



SPA AND WELLNESS CENTRES

Guests undergo profound rejuvenation and achieve a meditative reset. Health and well-being are the top priorities for travellers. Best suited for wellness centres, hotel spas and airport lounges.



MEDICAL AND TREATMENT FACILITY

Vibroacoustic Therapy can easily be integrated into rehabilitation centres, medical spas, and traditional psychology practices. It helps rebalance the nervous system, reduce stress, anxiety, and pain, and reduces inflammation in the body.

SPECIFICATIONS

FEATURES

Design	Zero Gravity positioning for a deeply relaxed, weightless feel
Foam Cushion	155 mm mixed layer of different densities for comfortable massage experience

TECHNICAL SPECIFICATIONS

Length	1915 mm / 75 in
Width	745 mm / 29 in
Height	772 mm / 30 in

ELECTRICAL PARAMETERS

Power in (Mains Power Supply)	100–240V, 50–60Hz, 3A
IPX Rating	IPX6
Operating Temperature	10–40°C
Shock Protection:	Class I

BASE & FRAME

Base	Powder-coated mild steel frame
Safe Working Load	250 kgs
Load Bearing capacity	180 kgs

SAFETY STANDARDS

General Safety Standards	IEC 60601-1
EMI/EMC	IEC 60601-2

UPHOLSTERY

Type	Phthalate-Free Upholstery
Thickness	1 mm
Tensile Strength	ISO 1421, 869.22 N/ 5 cm
Colour Options	Cream, Beige, Mouse Grey, Brown, Charcoal Grey
Weight	726.7 g / sqmt
Abrasion	ISO 5470-2 = 2,00,000 cycles / min

WARRANTY

Motors	2 years
Frame	5 years



Four Tactile Sound Transducers

The Harmony Vibroacoustic Lounger features four premium tactile sound transducers embedded directly within the table, converting sound into full-body vibration. Rather than simply hearing the wellness programs, guests feel every frequency resonating through the body, deepening the therapeutic effect of each session.

Six Pre-Loaded Wellness Programs



The Harmony comes with six expertly curated programs: Relaxation, Stress Relief, Sleep, Focus and Mental Clarity, Tibetan Gong and Singing Bowl, and a Voice Guided Demo. Each program is precisely tuned to target specific brainwave states, giving every guest a deeply personalised wellness journey.

Zero Gravity Wave Shape Ergonomic Design


The Harmony's distinctive wave-shaped top is engineered to cradle the body in a naturally supported, fully reclined position. This allows guests to surrender completely to the experience and enabling vibrations to travel through the body with greater depth and effect.

Phthalate-Free Upholstery


The Harmony is available in Phthalate-free Upholstery, designed for professional spa and wellness environments, offering a safer surface for both clients and therapists. It resists acetone and disinfectants, withstands frequent cleaning, and ensures durable, fire and abrasion-resistant performance.

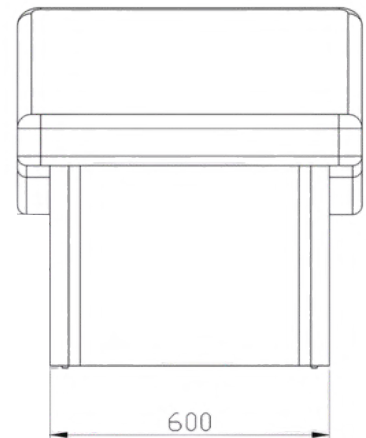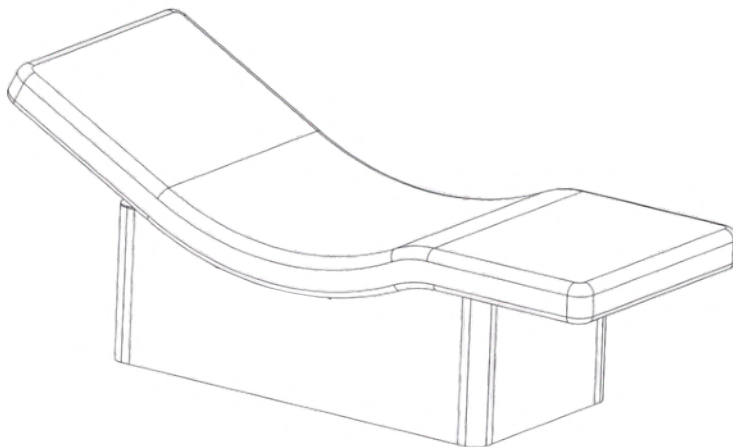

FRONT VIEW

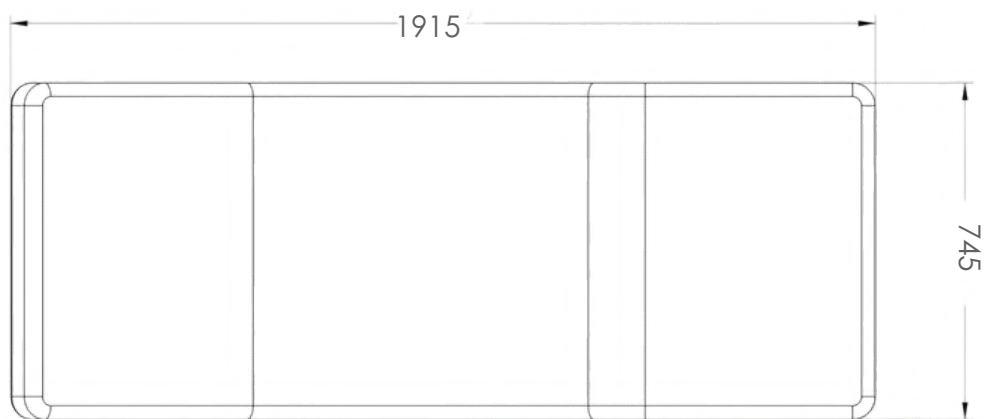

TOP VIEW

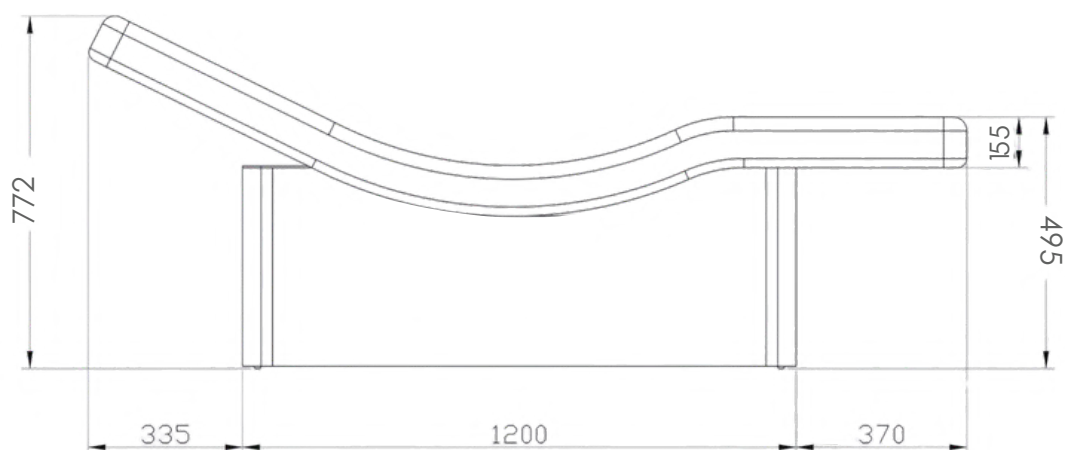

SIDE VIEW

ALL DIMENSIONS ARE IN MM

01.



Relaxation - 22 Minutes

- Alpha Brain Waves

Designed to promote a deep state of rest and calm.
02.



Stress Relief - 24 Minutes

- Theta Brain Waves

Intended to support a quieter, more settled state of mind.
03.



Tibetan Gong & Singing Bowl - 22 Minutes

- Alpha Brain Waves

An ancient-inspired soundscape to encourage stillness and meditative calm.
04.



Focus & Mental Clarity - 25 Minutes

- Gamma Brain Waves

Designed to support concentration and a sense of mental presence.
05.



Sleep - 23 Minutes

- Delta Brain Waves


A slow, deep program intended to encourage the mind toward restorative rest.
06.



Voice-Guided Demo - 10 Minutes


- Alpha Brain Waves

An introductory session ideal for first-time guests.

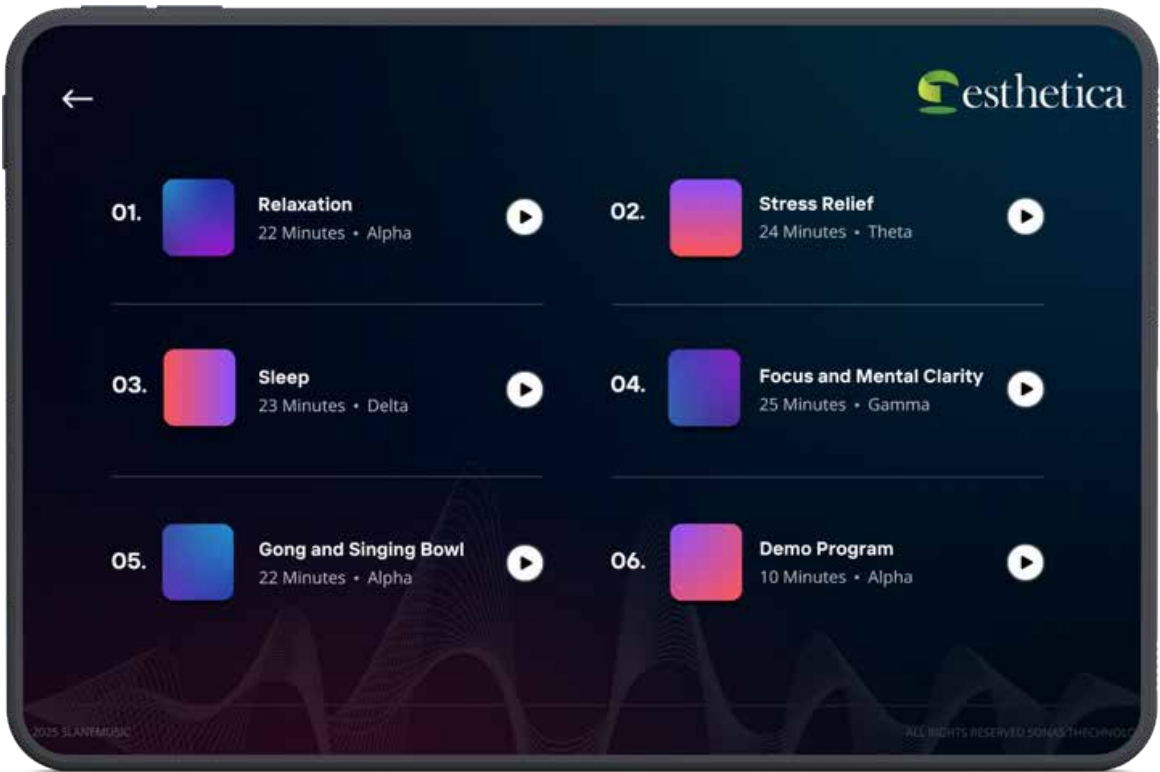



- Connect the tablet to the ESI Gigaport using a USB-C to USB-C cable.
- Plug in the tablet's Adapter and connected to power before turning it on. The tablet does not work on battery.

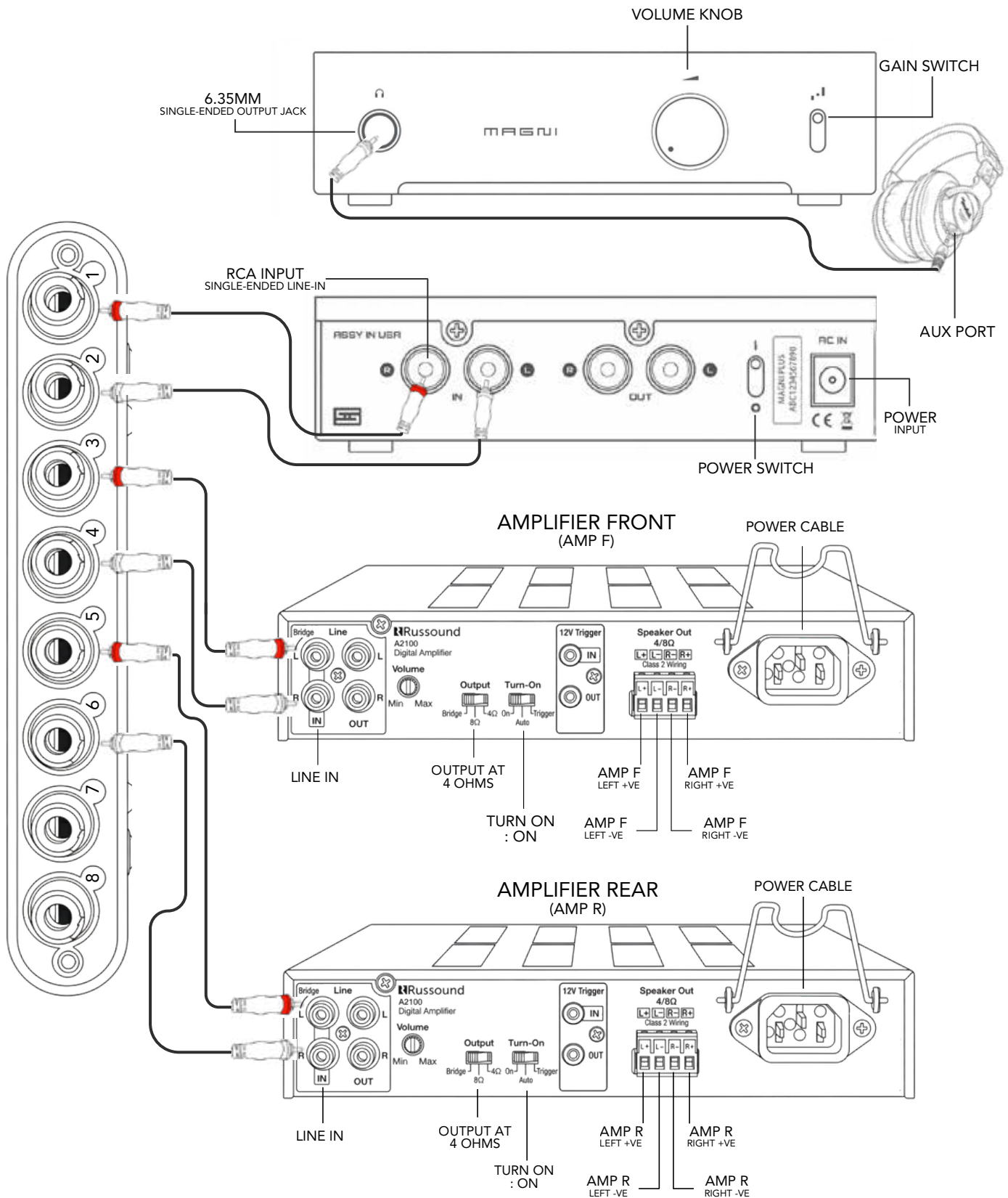


MAIN AMPLIFIER CONNECTIONS (RCA CABLES)

Using the RCA cables, connect the amplifiers to the ESI Gigaport as follows:


ESI Channel 1 & 2 → Line In of Headphone Amplifier

ESI Channel 3 & 4 → Line In of Front Amplifier

ESI Channel 5 & 6 → Line In of Rear Amplifier

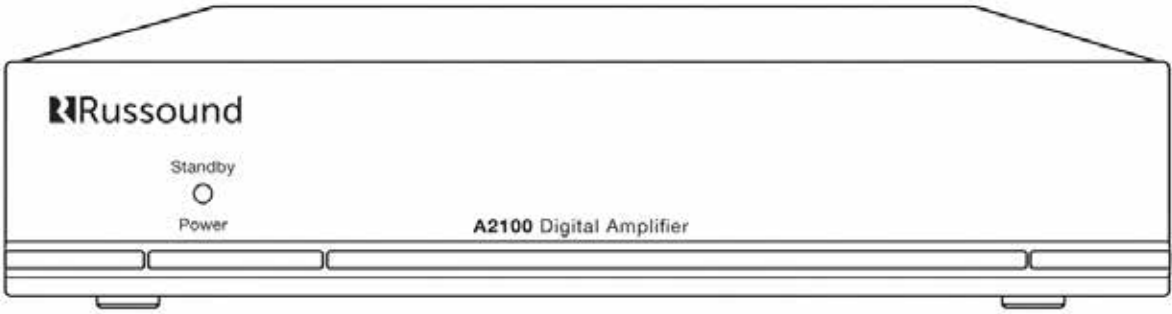

FRONT VIEW OF AMPLIFIER

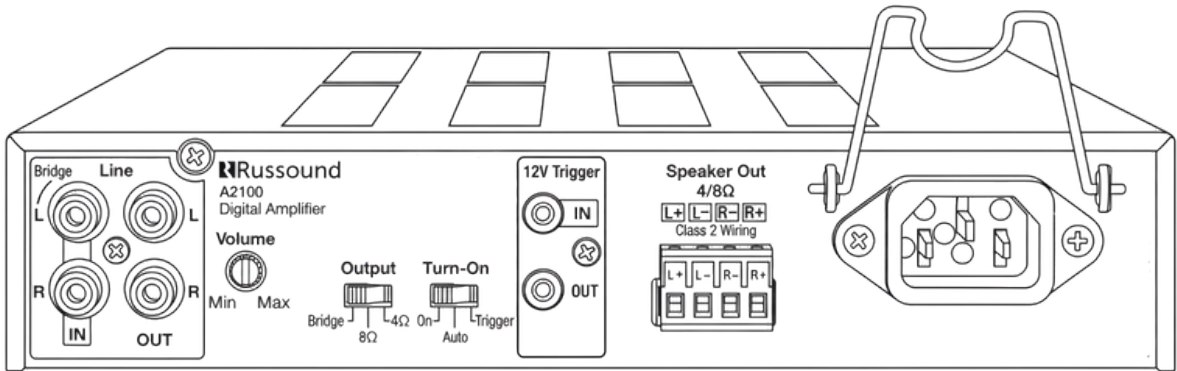

REAR VIEW OF AMPLIFIER

SPEAKER CONNECTIONS

Connect the speaker wires from the vibroacoustic bed to the Speaker Out of both main amplifiers:

- The amplifier connected to Channels 3 & 4 of the ESI Gigaport Ex (Front Amplifier) must be wired to the speaker leads labelled AMP F.
- The amplifier connected to Channels 5 & 6 of the ESI Gigaport Ex (Rear Amplifier) must be wired to the speaker leads labelled AMP R.

ESI GIGAPORT EX CHANNELS	AMPLIFIER	SPEAKER WIRE LABEL
3 & 4	Front (AMP F)	AMP F - LEFT+VE, LEFT-VE, RIGHT-VE, RIGHT+VE
5 & 6	Rear (AMP R)	AMP R - LEFT+VE, LEFT-VE, RIGHT-VE, RIGHT+VE

- Set the output of both main amplifiers to 4 Ohms.
- Plug in the power cables of both main amplifiers.
- Power on the amplifier using the front panel switch.

HOW TO CARE FOR YOUR UPHOLSTERY

In an effort to ensure your product quality and help you maintain the life of your table, we have tested several preventive steps in order to advise you regarding how to protect the upholstery on your table and accessories.

- Avoid storing your table in extreme cold temperature as excessive coldness will cause cracking.
- Avoid storing your table in extreme hot temperature as excessive heat will cause the upholstery to stretch and pucker. Be sure to keep your table away from heaters, electrical appliances and direct sunlight.
- Protect your table and accessories with sheets and fleece.
- Clean your upholstery daily.
- We recommend cleaning your upholstery with a mild detergent to protect it against oil stains and dirt, which will cause cracking and permanent stains.
- We recommend mild mixture of anti-bacterial dish soap with water in spray bottle. Spray the table with the solution and wipe down the table with a clean damp towel.
- Use both ends of the table. Occasionally, switch the face rest from one end of the table to the other to avoid overuse.
- Protect your table from general abrasions, cuts and punctures.

CAUTION: Prolonged exposure to extreme temperatures may cause damage to vinyl.



Phthalate-Free



Fire Retardant



Oil Resistant



Alcohol Resistant



Ultra Soft



Anti Bacterial



Eco - Friendly

Wood Finish Options - American Oak



American Oak



White




White Wash Finish




American Oak (Gloss Finish)



Teak

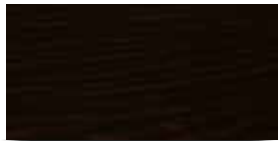

Wenge

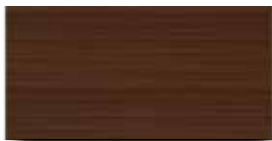


Dark Walnut




Black


Phthalate-Free Upholstery Options

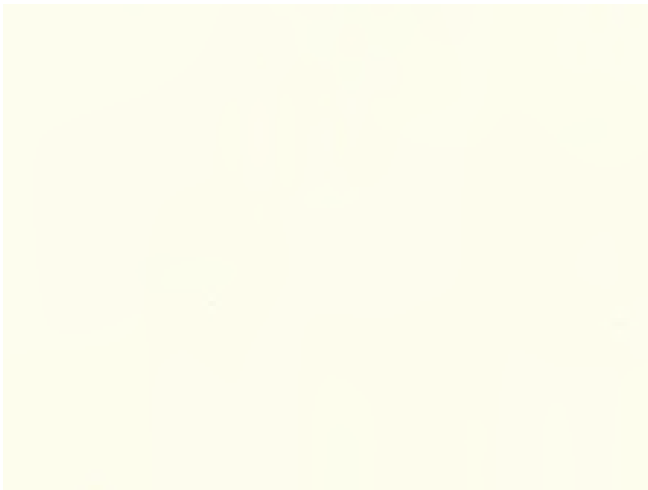


Cream

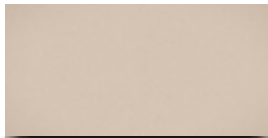


Beige

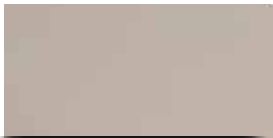

Mouse Grey

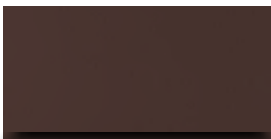

Brown

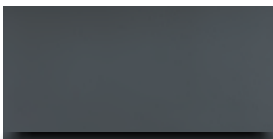

Charcoal Grey



Phthalate-Free



Fire Retardant



Oil Resistant



Alcohol Resistant



Ultra Soft



Anti Bacterial



Eco - Friendly